



# ARYS™ pro

Getting better. Everyday.



## What is ARYST™ pro?

ARYST™ pro is a professional medical device system for the documentation of the treatment of patient's arm and hand function due to hemiparesis after a stroke or traumatic brain injury. This system allows patients to automatically record even minimal arm movements. The ARYST™ pro system is based on the recording of all arm movements by means of two trackers, one of which is worn on the wrist of the non-paretic arm and one on the paretic arm and an innovative evaluation software running on an iPad tablet.

- ▶ Allows patients to follow daily progress thanks to a 4-LED status bar.
- ▶ Monitors arm activities and reminds the patient discreetly when paretic arm is neglected.

## Follow your patient with ARYST™ pro beyond the therapy sessions.



Ratio paretic vs. non-Paretic arm usage.

Arm activities:  
 • Blue: paretic  
 • Grey: non-paretic

Total of reminders for inactivity and step counter.


## ARYS™ pro – Everything on the tip of your finger!



- 1 Therapist profile
- 2 Trackers setting
- 3 Navigation bar
- 4 Data from last sync
- 5 Tracker battery status
- 6 Patient's details
- 7 Ratio paretic/non paretic
- 8 Improvement in paretic arm usage
- 9 Goals set
- 10 Latest note and assessments



## Reward your patient's efforts!

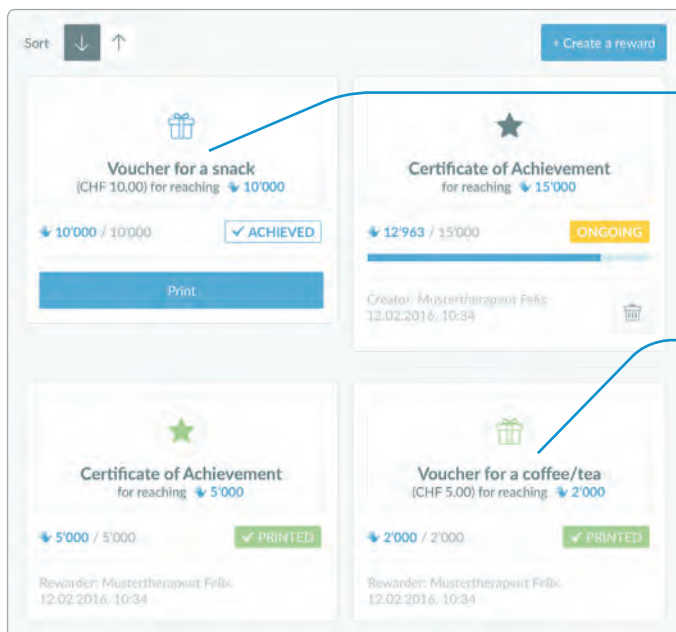

When the patient reaches the individual therapeutic goal set by the therapist, rewards can be given in the form of vouchers.

Look back and celebrate therapy goals reached – and hand over a personalized certificate of achievement. Success made visible has tremendous power! It cheers the patient to go on further.

- ▶ We help you and your patients to achieve the set goals while staying motivated.



Get more information at  
[www.ybandtherapy.com](http://www.ybandtherapy.com)

ARYS™ pro has been developed in close collaboration with leading neurologists and therapists in the field of neurorehabilitation. ARYS™ pro has been designed and manufactured according to the applicable standards. The tracker is protected against splash water.

### Manufacturer

**yband therapy AG**  
Laufenstrasse 66  
4053 Basel, Switzerland  
+41 61 515 64 80  
info@ybandtherapy.com

[www.ybandtherapy.com](http://www.ybandtherapy.com)



 yband  
therapy  
We empower rehabilitation.